


Model no. 3S-DAC90



Email: <u>info@3smedical.co.kr</u> Homepage: <u>www.3smedical.com</u>

Super Dried Air Compressor

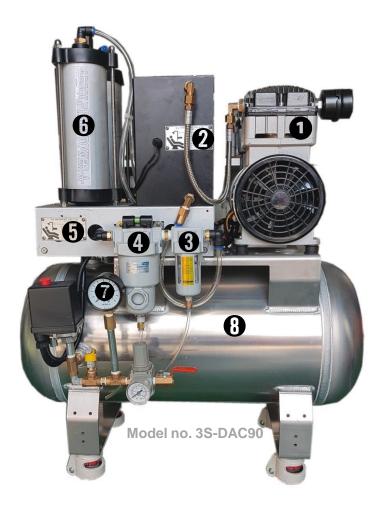

Description

1	Oil-less Air Compressor
2	Fin Coil type Air Cooler
3	Automatic Water Separator
4	Auto Drain-Bacteria Filter
5	Artificial Intelligent clean-air generator
6	Dual operation Air Drier
7	Pressure controller
8	Aluminum storage Tank

Air Compression



Air Cooling



Automatic Water separation

Our Process

Aluminum Storage tank



Dual operation
Air drying



Removed Germ



WHAT IS DIFFERENT?



Membrane Air Dried Compressor



- High Cost Membrane FilterShort-term Replacement of Filter
- Complicated Maintenance

■ High Cost Membrane Dryer


- Reasonable Cost Adsorption Dryer
- Use Semi- permanently Drying Media
- Easy Maintenance



CLINIC SMELLS STINKY?



- **COMPETITOR** RING BLOWER & WATERING PUMP METHOD

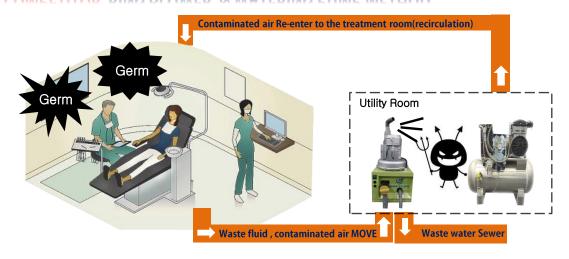

-3S MEDICAL VACUUM COMPRESSOR METHOD



In many dentists, fatal environment is given to patient, doctor, and medical staff everyday (30 times of general bacteria, 5 times of lead, 3 times of fluorine, 11times of ammonia nitrogen occur)

3S Medical products completely separate the absorbed saliva and waste solution with specialized dual separator, even separated odor is discharged to the drain with waste water. In conclusion, we improved hygiene condition for all the people in the dentist.